

Innovative Technology
For Your Medical Aesthetic Needs

COOLSHAPING²

Extreme Cooling Lipolysis device

COOLSHAPING²

SPECIFICATIONS

Energy type	Cooling lipolysis
Cooling range	-10~-5°C
Vacuum pressure	260mmHg

Handpiece	Large, Medium, Small, Mini
Input power	AC 100~240V, 50 / 60 Hz
Dimension (WxDxH)	498 x 774 x 1,119 (mm)


 **Eunsung**

Head office : 120, Gieopdosi-ro, JiJeong-myeon, Wonju-si, Gangwon-do, Korea. Tel : +82-33-760-1700(Rep.) Fax : +82-33-746-9716
Seoul office : #1701-1702, 145, Gasan digital 1-ro, Geumcheon-gu, Seoul, Korea. Tel : +82-2-514-9713(Rep.) Fax : +82-2-514-9716
E-mail : es@esglobal.co.kr / www.esglobal.co.kr



Treatment by cooling lipolysis technology

Coolshaping2 is a cooling lipolysis device. Combination of Cryo and Vacuum technology creates a safe and secure body-sliming effect. Compared to the existing device, the treatment area has been enlarged and the treatment time has been greatly reduced. Various types of 4 handpieces allow the treatment of local areas such as face.

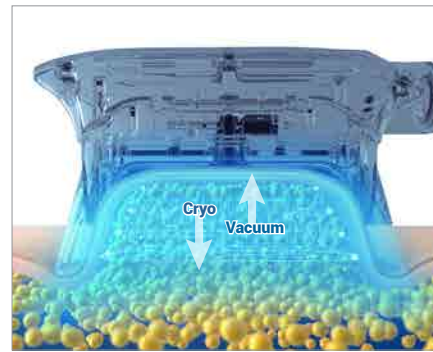


Mini handpiece
*70 x 139 x 200 (mm)

Small handpiece
*89 x 139 x 119 (mm)

Medium handpiece
*93 x 188.7 x 96 (mm)

Large handpiece
*115 x 225.5 x 120 (mm)

Treatment Area

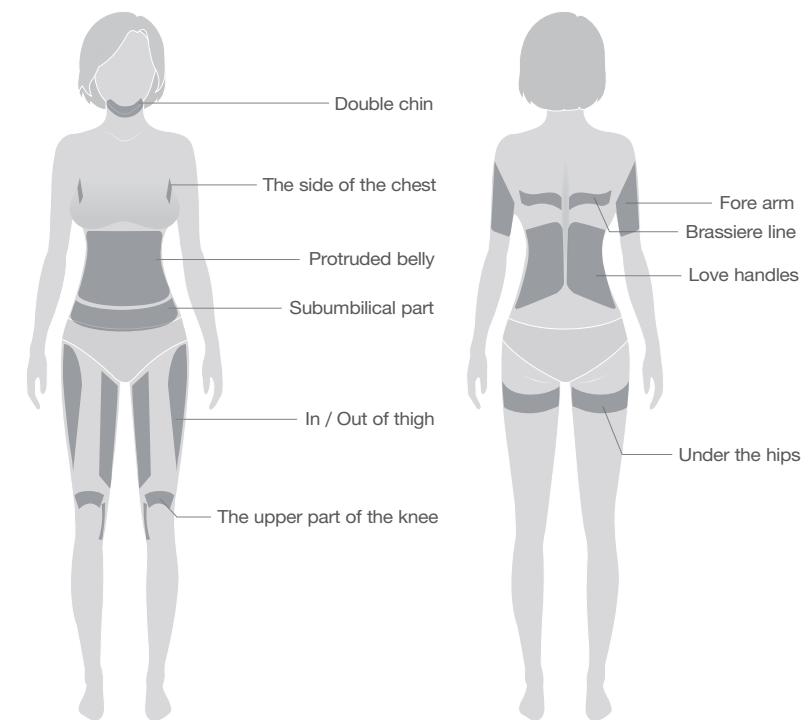

Coolshaping's Revolutionary Change!

COOLSHAPING²

- ✓ Now, 360° All Area Cooling!
- ✓ Simultaneous operation of 4 handpieces!
- ✓ 4 Types of Handpiece for various treatment areas!
- ✓ Enlarged treatment area!
- ✓ Greatly reduced treatment time!

360°

Treatment Mechanism



Start Cryo-cooling

Apply coolshaping applicator at the target area, start cryo-cooling to reach target temperature.



Immediate Cryo-reaction

The fat cells in the target area begin responding to cryo-crystallization.



30 Days Metabolism process

Over the following weeks and months, fat cells are disrupted and eliminated gradually by metabolism.



90 Days Eliminated cells

The fat cells are eliminated.

Before & After



Before



After

Abdomen F / 32
3 treatments



Before



After

Abdomen M / 40
2 treatments



Before



After

Thigh F / 26
3 treatments

Advantages



Stable Cryo energy



Various applicators



Shot Downtime



Simultaneous operation



User convenient interface

Accessories



◀Cushion



◀Cool pad



◀Liner