

Innovative Technology
For Your Medical Aesthetic Needs

FOCUS_{DUAL}

RF & HIFU
Multi-System

Skin tightening
Acne scars treatment
Anti- wrinkle
Skin lifting
Skin brightening

No Discomfort
No Side Effect

FOCUS_{DUAL}

SPECIFICATION

| | |
|-------------------------|--|
| Energy Type | Radiofrequency & HIFU |
| Handpieces | RF Handpiece & HIFU Handpiece |
| Catridge | 25 pin SM4-4,5 / DD7-3,0 / SD7-1,5 SD7-1,5N / SC4-13 |
| Input Power | AC 100~240V, 50/60Hz |
| Dimension(W x D x H mm) | 504 x 438 x 329 |



 Eunsung

Head office : 120, Gieopdosi-ro, Jijeong-myeon, Wonju-si, Gangwon-do, Korea. Tel : +82-33-760-1700(Rep.) Fax : +82-33-746-9716
Seoul office : #1701-1702, 145, Gasan digital 1-ro, Geumcheon-gu, Seoul, Korea. Tel : +82-2-514-9713(Rep.) Fax : +82-2-514-9716
E-mail : es@esglobal.co.kr / www.esglobal.co.kr

RF & HIFU 2-in-1 System

FOCUS DUAL

Combination of Fractional Micro Needle RF Technology & HIFU
The best solutions for divers skin complex is Focus Dual!

There are many diverse skin problems. With only one treatment, skin complex can not perfectly managed. Focus dual is the combination device of fractional micro needle RF & High intensity focused ultrasound technology. The device can be used from acne scar treatment to lifting effects. The two handpieces on Focus Dual are both non-invasive for the epidermis. Non-invasive treatment is safer way for both operator and patient. Especially, the delicate setting adjustment will give the operators high satisfaction with convenient use and significant clinical effects.



Fractional Micro Needle RF & HIFU




| | Fractional micro needle RF | HIFU |
|-----------------------|---|---|
| Precision | <ul style="list-style-type: none"> Depth : 0.5-3.5mm RF on time is able to control | <ul style="list-style-type: none"> RF level : 1-10 level Suction : 0-3 level |
| Advantages | <ul style="list-style-type: none"> Able to select either foot switch or handpiece switch. Safer treatment by applying vacuum function. Hygienic treatment for using disposable multi needle. | <ul style="list-style-type: none"> Easy cartridge replacement. Comfortable 3 buttons design 1 cartridge=20,000 shots |
| Comfortability | <ul style="list-style-type: none"> 10.2 inch touchscreen User friendly UI settings Ergonomically designed handpiece able for long-time use | |

Fractional Micro Needle RF

When micro-needle is injected into the skin, on the same time, Bi-polar RF will be also produced. The RF energy will be directly affected to the dermis and help regeneration of collagen and elastin. Also, needling can improve acne scars and reduce pores. By adjusting proper needling depth, it can solve diverse skin problems.

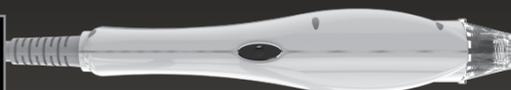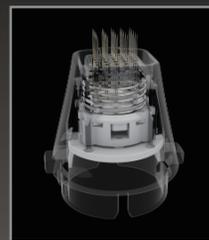



- 33G 25pin (5x5) ,non-isolated needle
- Non-isolated needle: Apply on the proper skin layer based on each target
- By suction(Vacuum) technology, able to do the treatment with proper depth to the skin
- RF energy, which able regular needling interval
 - By applying Rf energy without overlapped needling, it can offer definite fractional effects

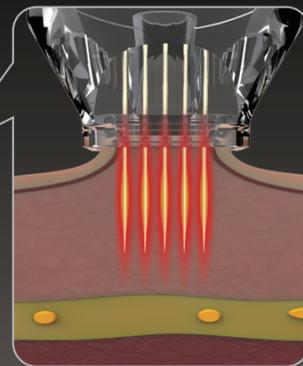

Effects Before & After



HIFU High Intensity Focused Ultrasound

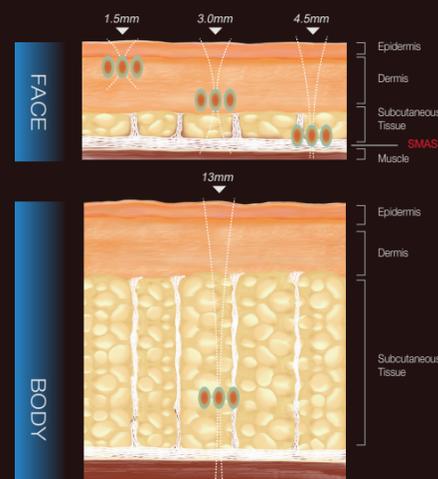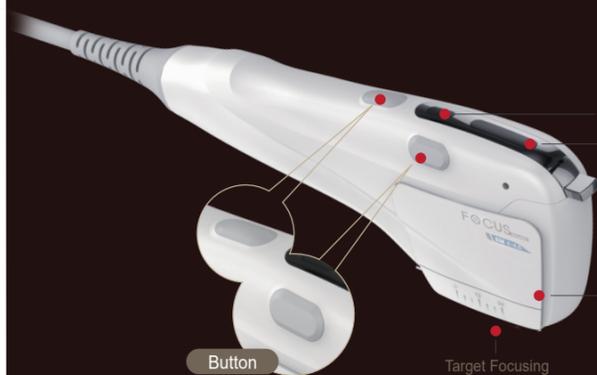

High intensity focused ultrasound technology is a non-invasive treatment that has no damage on the epidermis. There are four cartridges for facial-use. Through the treatment, the elasticity of the skin is improved and has lifting effects by concentrating the ultrasound energy on each layer of the skin, such as upper dermis, lower dermis or SMPS. Moreover, there is also a cartridge for the body, with 13mm depth (subcutaneous layer) to reduce the fat effectively.

